



AirFit™ F20
AirTouch™ F20



Nine times quieter and just as effective

Introducing QuietAir™, a revolutionary diffuser vent that makes ResMed AirFit™ F20 and AirTouch™ F20 quieter than ever before.¹

Compared with our standard vent, QuietAir delivers the same level of performance while minimizing noise and disruptive drafts.

[ResMed.com/QuietAir](https://www.ResMed.com/QuietAir)




9x quieter
than before*

70% gentler
exhaled air flow†

Performance factor	QuietAir vent vs. standard vent
CO ₂ released	Equal
Vent flow	Equal
Resistance to flow	Equal
Quietness	9x quieter
Airflow	70% gentler airflow

QuietAir disperses the same amount of air per breath as the existing standard vent, but it does so using smaller, multidirectional openings, which spread the exhaled airflow over a wider area. The result is airflow so soft that patients might not feel it at all.

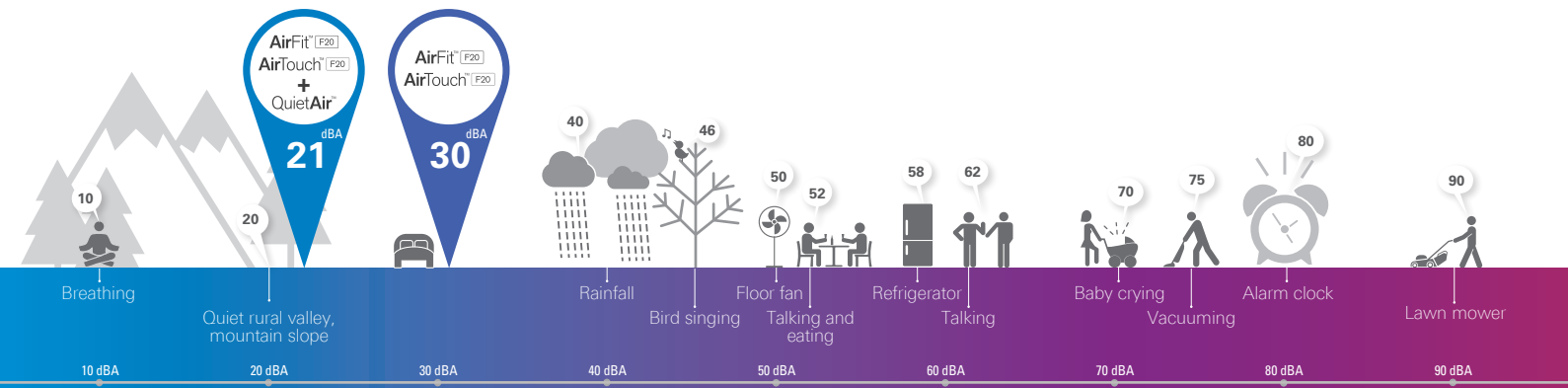

Reference: Berger EH, Neitzel R and Kladden CA. 3MTM Noise Navigator™ Sound Level Database Version 1.8. 22 August 2017. Available at: <http://multimedia.3m.com/mws/media/888530/noise-navigator-sound-level-hearing-protection-database.pdf>

1 ResMed, 2017. [Studies and claims for AirFit F20, ID 1019565]. Unpublished internal document.

* Reduction in sound power level of 9 dBA compared with the standard F20 elbow.†

† Measured as the reduction in peak velocity when compared to the standard F20 elbow.†